



UHV Cryostat - Ultra Low Vibrations

The CS204*F-DMX-20B opens up a wide range of possibilities never before possible with a closed cycle cryocooler. The ARS manufactured DMX-20B interface, like our “B” series cold heads is a True UHV System (10^{-11} Torr) where all of the rubber o-ring seals have been replaced with welded joints and metal seals. A CF flange is The DE204*F-DMX-20B uses a Helium Exchange Gas to decouple the sample from the cold tip of the DE204 Cryocooler. This prevents almost all vibration from being transmitted to the sample. Sample vibration levels of 3-5nm have been demonstrated by users. Due to the exchange gas being less conductive, the base temperature will increase by 1-2K and the cooling capacity is roughly cut in half at varying temperatures.

Applications

- Low Vibration UHV applications
- Microscopy Applications
- Surface Science
- Nanomaterials

Features


- Ultra Low Vibrations (3-5 nm)
- True UHV (10^{-11} Torr)
- Bakeable to 200C with cold head removed
- Open Sample Space
- Cold Tip Down Orientation
- Fully customizable

Typical Configuration

- Cold head (DE-204AF)
- Compressor (ARS-4HW)
- 2 Helium Hoses
- 6” Rotatable CF flange
- Nickel Plated OFHC copper radiation shield terminating 0.125” short of the cold tip
- 2 High purity quartz windows
- Instrumentation for temperature measurement and control:
 - 10 pin hermetic feed through
 - 50 ohm thermfoil heater
 - Silicon diode sensor curve matched to ($\pm 0.5K$) for control
 - Calibrated silicon diode sensor (± 12 mk) with 4 in. free length for accurate sample measurement.
- Wiring for electrical experiments:
 - 10 pin hermetic feed through
 - 4 copper wires
- Sample holder for optical and electrical experiments
- Temperature Controller

Options and Upgrades

- 4.5” CF flange available but uses non standard ID and may not fit existing chambers
- 8” CF flange available
- 4K Coldhead (0.2W @ 4.2K)
- 5.5K Coldhead (3W @ 10K)
- 450K High Temperature Interface (Not required with the DMX-20 interface, but the 450K interface is helpful for high temperature performance)
- 800K High Temperature Interface
- Turbo upgrade for faster cooldown times
- Custom temperature sensor configuration (please contact our sales staff)
- Custom wiring configurations (please contact our sales staff)
- Sample holder upgrades (custom sample holders available)

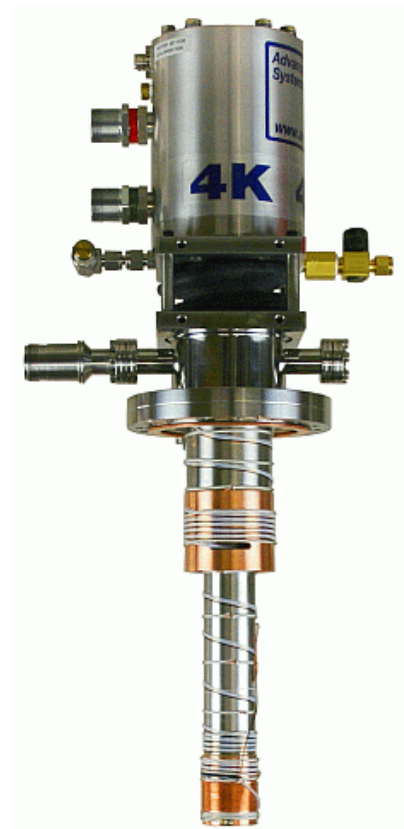

The above picture shows a cryocooler with a vacuum shroud, radiation shield, and sample holder installed.



The above picture shows the compressor for the system. What is not shown is a required vacuum pump, Helium Hoses (included with system), temperature controller and Helium Gas bottle with regulator for the exchange gas.



UHV Cryostat - Ultra Low Vibrations

Cooling Technology

DE-204	Closed Cycle Cryocooler
Refrigeration Type	Pneumatically Driven GM Cycle
Liquid Cryogen Usage	None, Cryogen Free

Temperature* DMX-20 adds 1-2K to base temperature

DE-204AF	< 10K - 350K
DE-204P4	< 4.2K - 350K
DE-202SF	< 5.5K - 350K
With 800K Interface	(Base Temp + 2K) - 700K
With 450K Interface	(Base Temp + 2K) - 450K
Stability	0.1K
*Based on bare cold head with a closed radiation shield, and no additional sources of experimental or parasitic heat load	

Sample Space

Diameter	Large Open Radiation Shield
Height	Large Open Radiation Shield
Sample Holder Attachment	1/4 - 28 screw
Sample Holder	www.arscryo.com/Products/SampleHolders.html

Chamber Interface

Flanges	CF, ISO
Size	4.5", 6", 8" 10" (4.5" CF flange non standard may not fit all)

Temperature Instrumentation and Control (Standard)

Heater	50 ohm Thermofoil Heater anchored to the coldtip
Control Sensor	Curve Matched Silicon Diode
Sample Sensor	Calibrated Silicon Diode with free length wires
Contact ARS for other options	

Instrumentation Access

Instrumentation Skirt	Welded Stainless Steel
Pump out Port	0
Instrumentation Ports	2
Instrumentation Wiring	Contact sales staff for options

Radiation Shield

Material	Nickel Plated OFHC Copper
Attachment	Threaded
Optical Access	Open End Radiation Shield terminates 0.125" short of cold tip (customer specified)


Cryostat Footprint

Overall Length	590 mm (23.24 in) 295 mm (11.63 in) standard flange to tip length
Motor Housing Diameter	114 mm (4.5 in)
Rotational Clearance	200 mm (8 in) with "G" Configuration

Cryocooler Model

		DE-204AF		DE-204A(T)F		DE-204PF		DE-204SF	
	Frequency	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz	60 Hz	50 Hz
Base Temperature		<9K	<9K	<9K	<9K	<5.5K	<5.5K	<4.2K	<4.2K
Cooling Capacity	4.2K	-	-	-	-	-	-	0.2W	0.16W
	10K	2W	1.6W	2.7W	2.2W	3.5W	2.8W	4W	3.2W
	20K	9W	7.2W	12W	9.6W	8W	6.4W	8W	6.4W
	77K	17W	14W	23W	18.4W	14W	11W	14W	11W
Radiation Shield Cooling Capacity		18W	14W	24W	19W	18W	14W	18W	14W
Cooldown Time	20K	30 min	36 min	25 min	30 min	40 min	48 min	40 min	48 min
	Base Temperature	60 min	72 min	50 min	60 min	80 min	102 min	90 min	108 min
Compressor Model		ARS-4HW		ARS-4HW		ARS-4HW		ARS-4HW	
Typical Maintenance Cycle		12,000 hours		8,000 hours		12,000 hours		12,000 hours	


DE202*F-DMX-20 Vibration Spectra

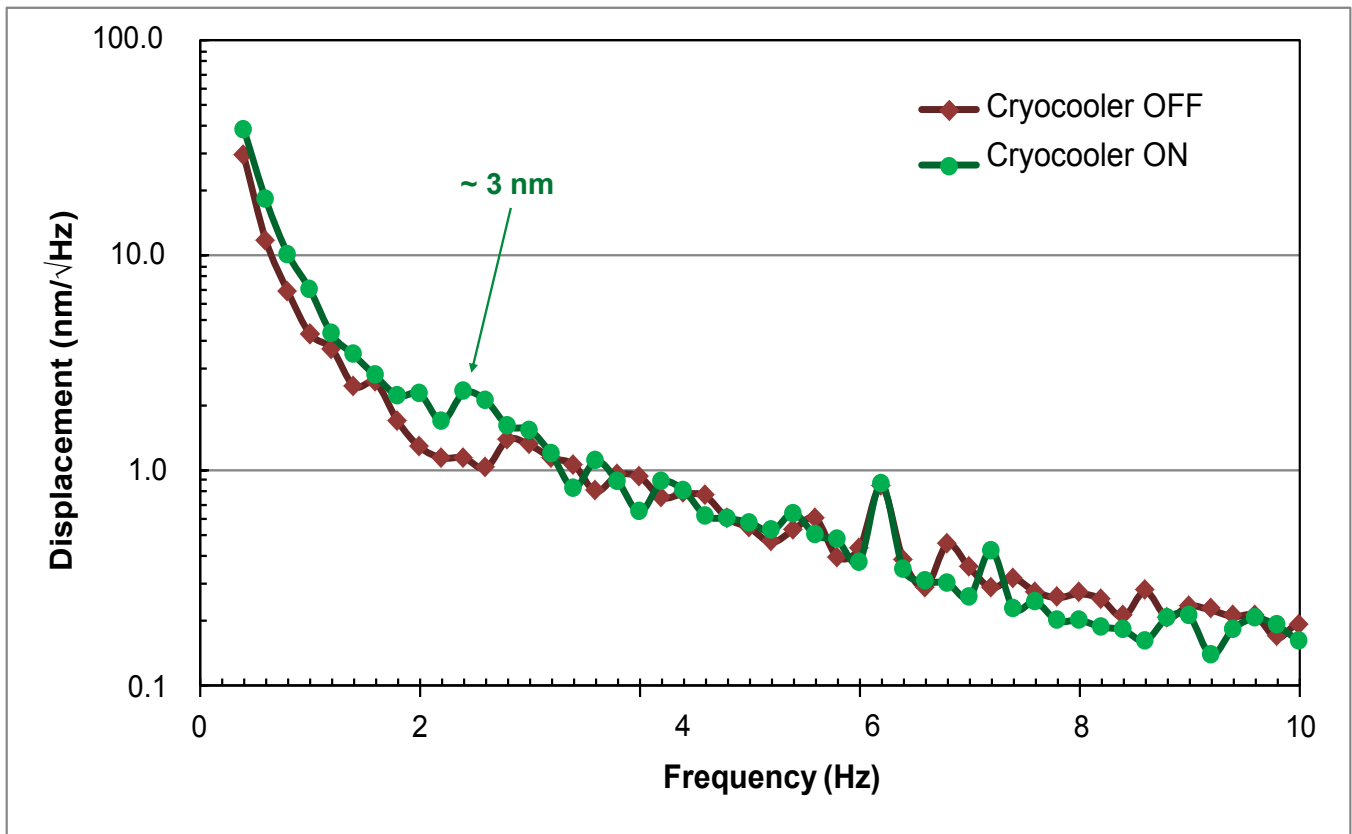

Understanding the DMX-20 Interface

The X-20 Interface uses a Helium Exchange Gas to decouple the sample space from the cold tip of the cryocooler. This prevents almost all vibration from being transferred to the sample space. Scientists have demonstrated vibration levels as low as 3-5 nm with the DE202*F-DMX-20 (as shown above)


1. The Cryocooler is supported from a Floor Stand
2. The soft rubber bellows minimize vibrations transmitted to the sample while keeping in the Helium Exchange Gas.
3. The cold tip has 10-30 micron vibrations (depending on CCR model) but no direct contact with the sample space.
4. Convective pockets of Helium Exchange Gas cools the sample space.
5. The sample is only in contact with the X-20 Interface
6. The X-20 Interface is mounted directly on a (user provided) Vibration Isolation Table.

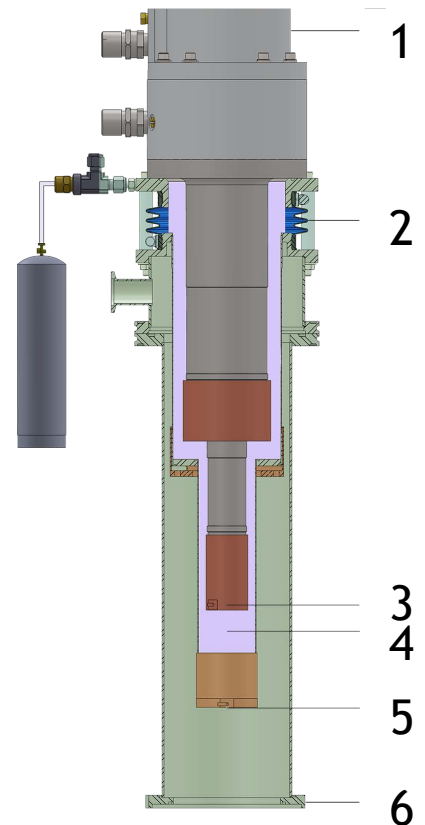



UHV Cryostat - Ultra Low Vibrations

DE204*F-DMX-20B Outline Drawing

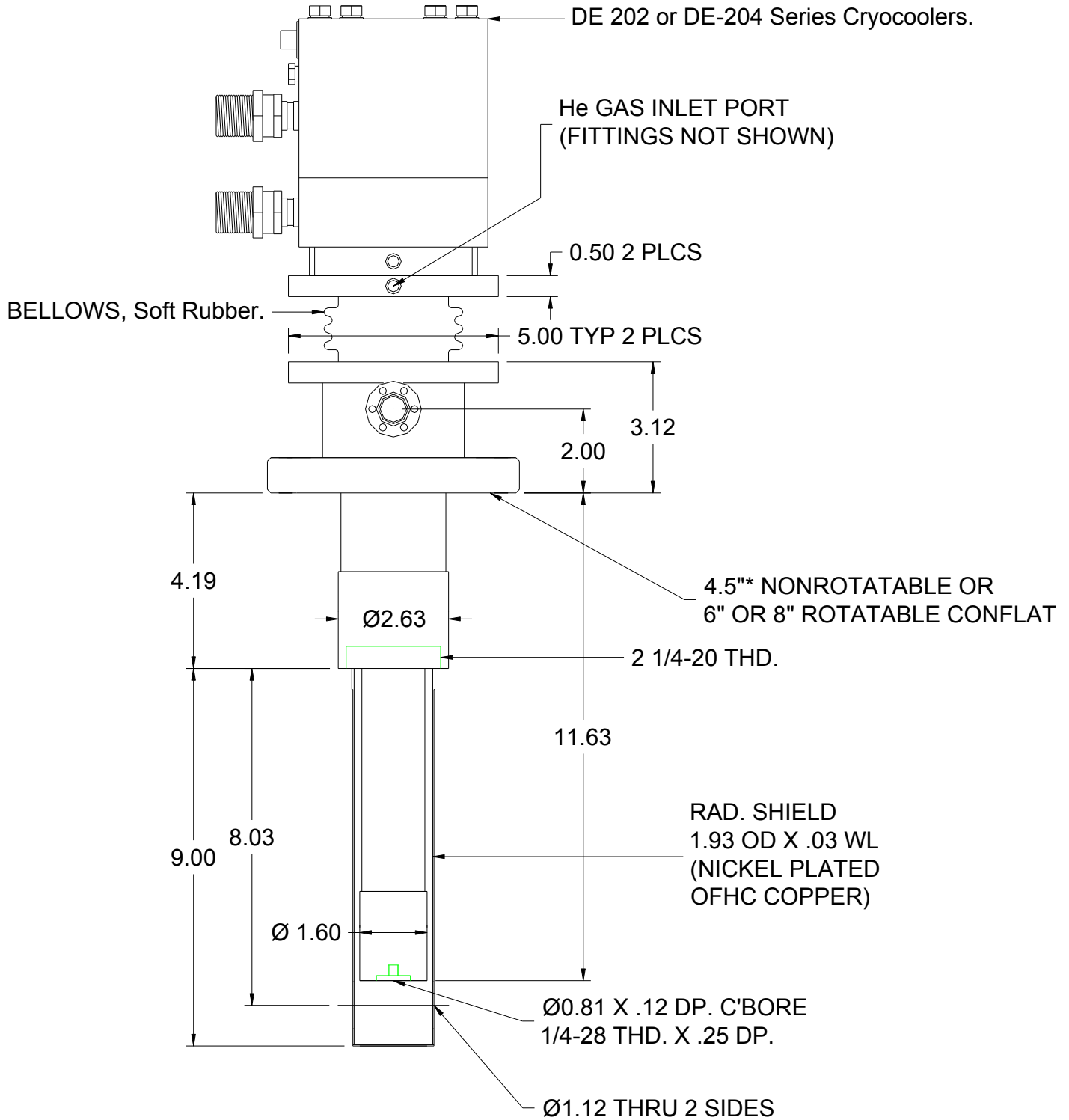




UHV Cryostat - Ultra Low Vibrations

Direct Mounting



The DMX-20 can be direct mounted on the Cryocooler. The vibrations at the sample will go up to 140 nm. It can be useful if the sample has to be translated in XYZ.

ARS-4HW Compressor

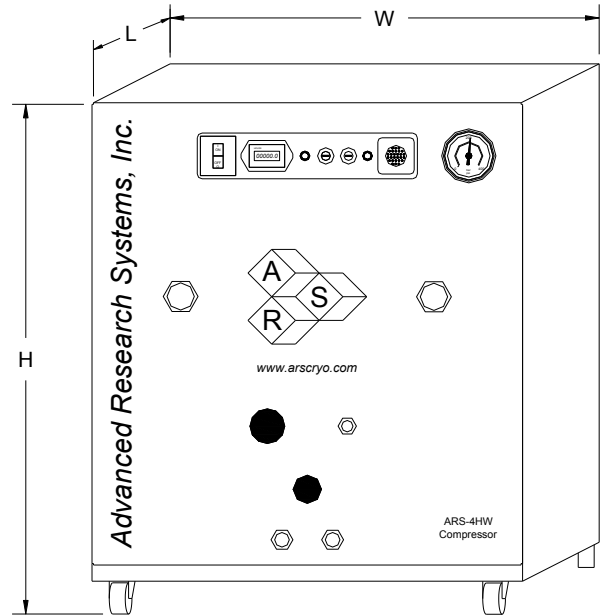


Compressor Model

ARS-4HW


	Frequency	60 Hz	50 Hz
Standard Voltage	Min	208 V	190 V
	Max	230 V	210 V
Transformer Options	10%	220 V, 230 V	
	15%	240 V	
Power Usage	Single Phase	3.6 kW	3.0 kW
Refrigerant Gas		99.999% Helium Gas, Pre-Charged	
Noise Level		60 dBA	
Ambient Temperature			
Cooling Water	Consumption	2.3 L / min (0.6 Gal. / min)	
	Temperature	10 - 35 C (50-95 F)	
	Connection	3/8 in. Swagelok Fitting	
Dimensions:	L	483 mm (19 in)	
	W	434 mm (17.1 in)	
	H	516 mm (20.3 in)	
Weight		72 kg (160 lbs)	
Typical Maintenance Cycle		12,000 hours	
Water Recirculation Option		CoolPac Compatible	

Floor Stand

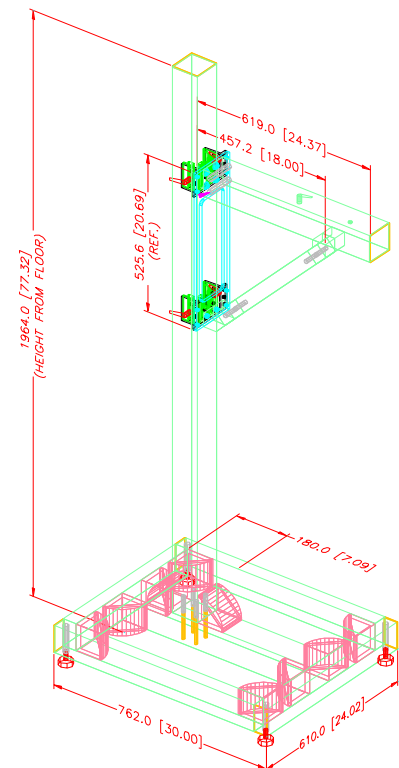

Mossbauer Spectra with DMX-20

Cryocooler Off

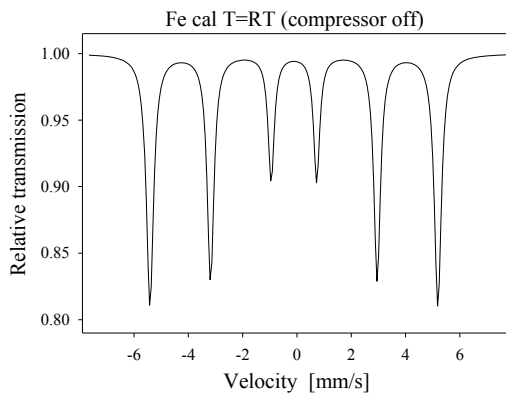

Calculated parameters:

WID=0.269 mm/s
W13=1.17

W23=1.08


ISO=0.000 mm/s

BHF=32.94T

Calibration spectrum:

Measurement with metallic iron foil (thickness 25 mm)

Cryocooler On

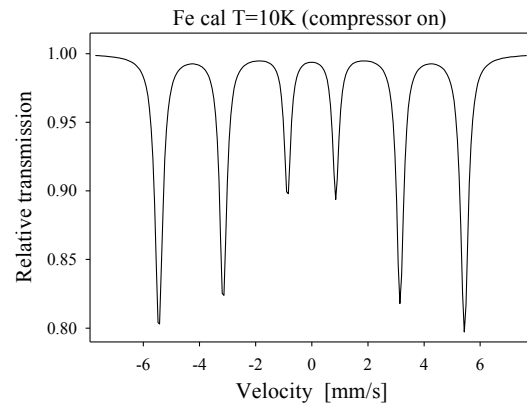

Calculated parameters:

WID=0.270 mm/s
W13=1.21

W23=1.11


ISO=+0.115 mm/s

BHF=33.81T

Calibration spectrum:

Measurement with metallic iron foil (thickness 25 mm)

Mossbauer Cryostat, DMX-20

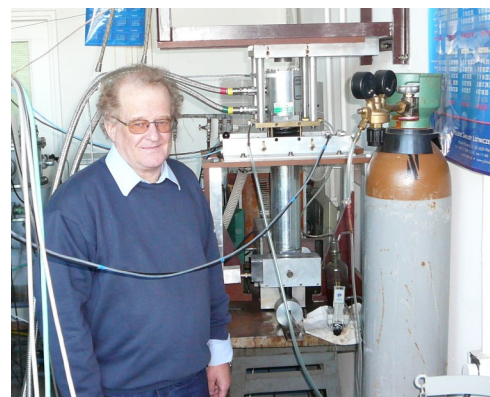

Prof. Dr. habil. Michal Kopcewicz

When the system is properly mounted to the wall and the table and properly adjusted then the effect of vibrations induced by the working compressor is very small (almost negligible).

The broadening (rather no broadening!) is fully reproducible and is well below 2%. Such broadening has practically no effect and in almost all cases can be ignored. However, compare the line widths calculated in the same way!

The Displex system is virtually vibration free!

Prof. Dr. habil. Michal Kopcewicz


Institute of Electronic Materials Technology

WARSAW, Poland



UHV Cryostat - Ultra Low Vibrations

SEM Setup

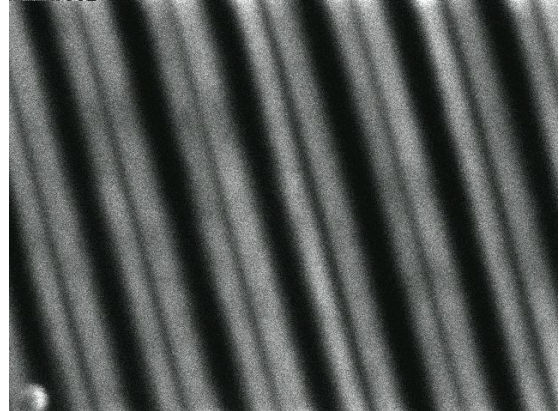

The picture shows the Displex and its interface. All items showing stainless steel in the photo are modifications to our JSM-5910 SEM.

Please also note the published papers on my website, which demonstrate CL results that were acquired with this cooling system.

The sample temperature, is between 40 K and 50 K.


Courtesy; Prof. Dan Rich, Ben Gurion University, Beer Sheva, Israel.

Image of an InGaN/GaN quantum well grown on a patterned GaN substrate. The patterning consists of stripes with a periodicity of 10 microns.

Monochromatic CL image (Mag. 2000) taken with a 390 nm detection wavelength.

The degradation of image quality from room temperature and 50K suggests that the vibration is less than 10nm.

FTIR Setup

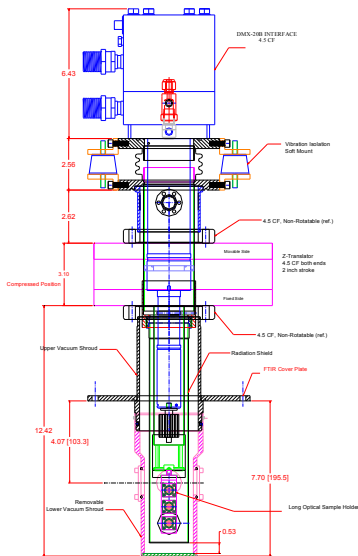


Low Vibration system for FTIR spectrometer. Sample holder with 3 samples can be translated in Z direction.

Low Vibration system for BOMEM, DA8 FTIR spectrometer.



UHV Cryostat - Ultra Low Vibrations


Magneto Electrical Experiments

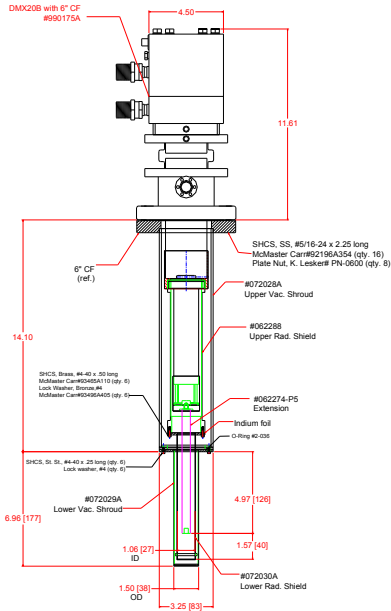

The vacuum shroud becomes narrow to permit sample insertion into a small magnet pole space.

The vacuum shroud becomes narrow to permit sample insertion into a small magnet pole space.

Magneto Optical Experiments (MOKE)

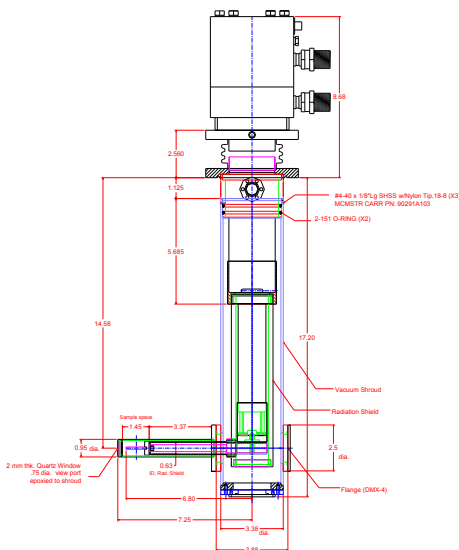

Low Vibration Side looking window can be placed in a MOKE, (Magneto Optical Kerr Effect). Sample can be located in any plane. The pole spacing can be as low as 1 inch.

Small diameter vacuum housing can be inserted into a narrow gap (high field). The optical window allows collection of light from the sample. The window can be very close to the sample for short focal length objective.